


**Step by Step HP Printer Installation Guide for Carl  
Zeiss IOL-Master 3.0 for Windows 98**


## HP Printer Installation Procedure

Unplug the Centronics printer cable from the “**DeskJet**” or other printer of the IOL-Master device. Connect the same cable to the Print2USB box.


From the IOL-Master main menu, select “**Options**”

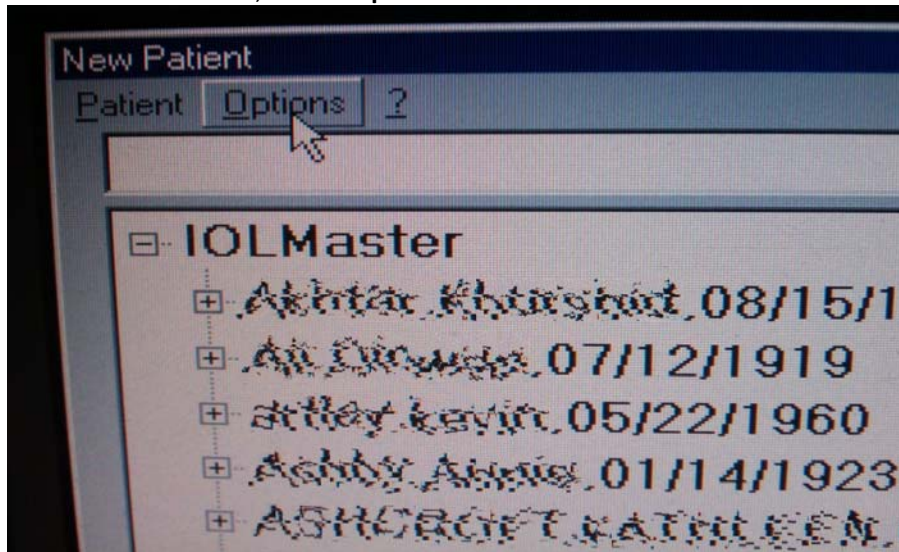


From the drop down menu select **Setup** -> “**Add Printer**”

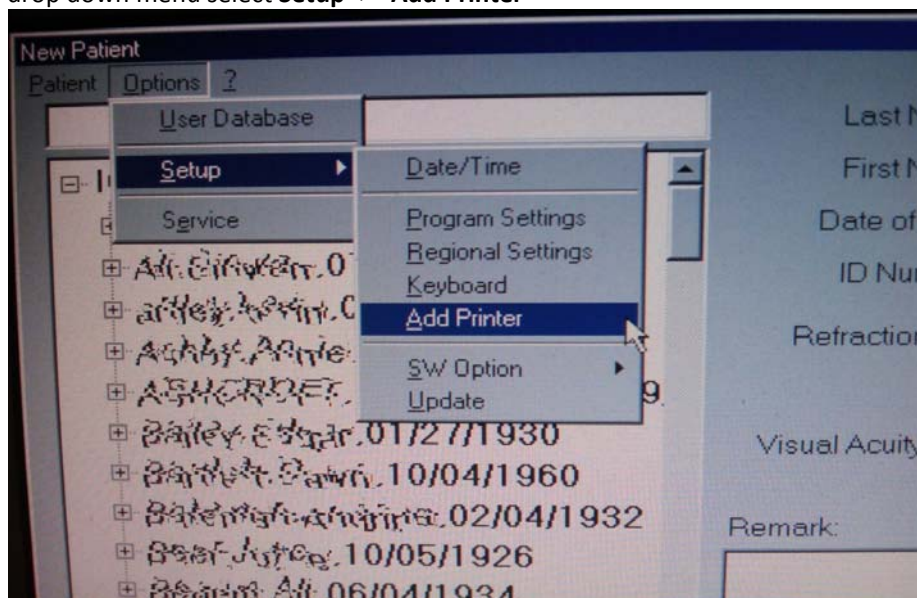


On the new window double click on “Add Printer” icon

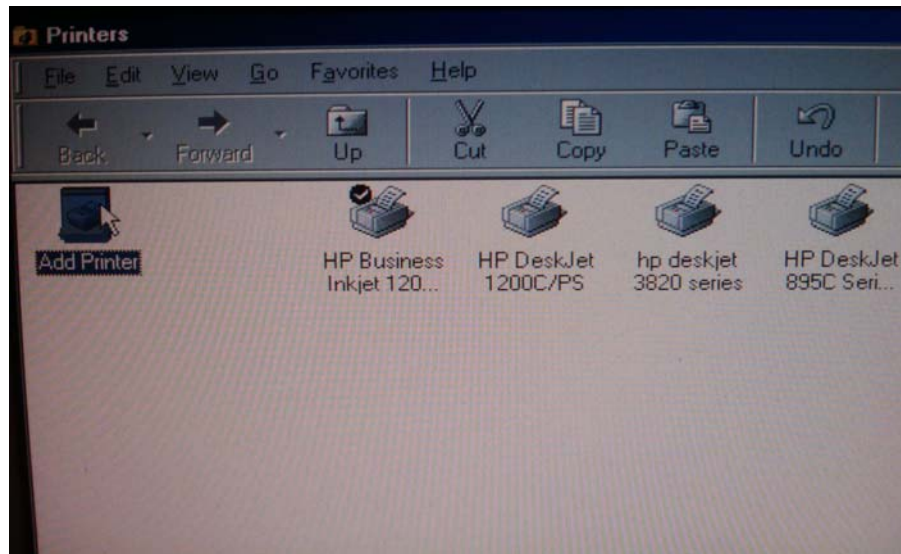


Click on “Next”

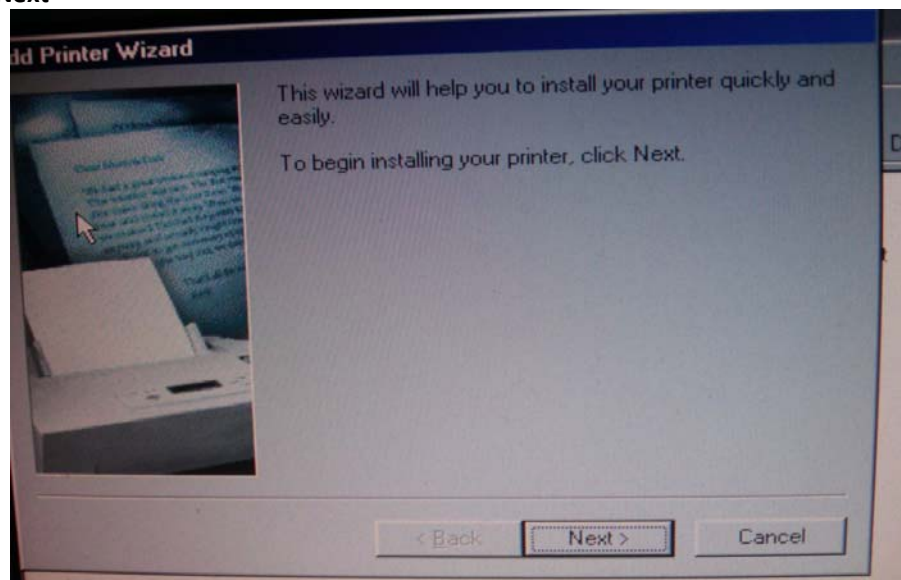


Select local printer option then click “Next”

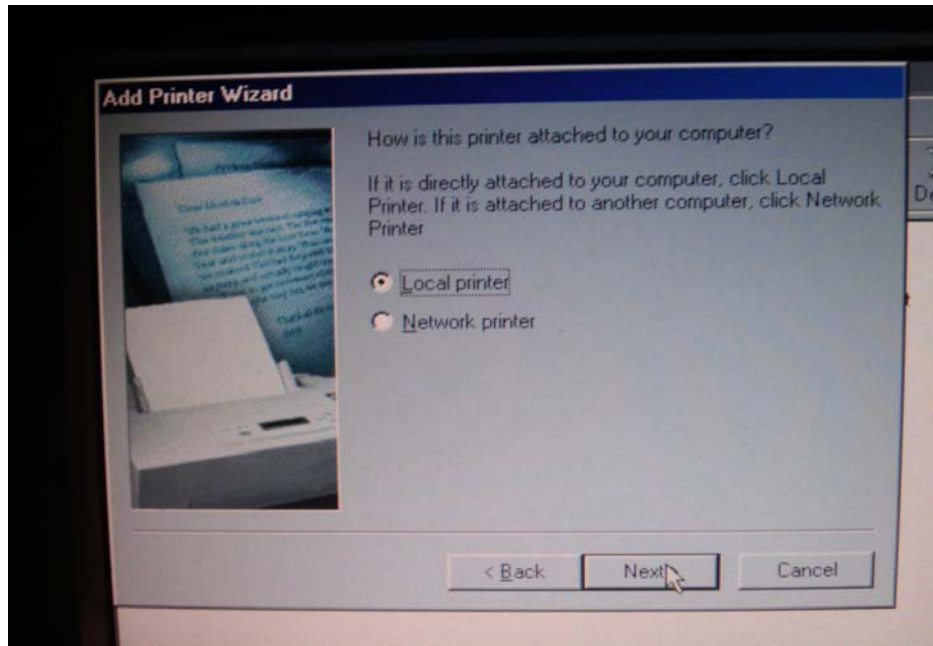


In the “add printer wizard”, Select “HP” under “Manufacturers” drop down box.

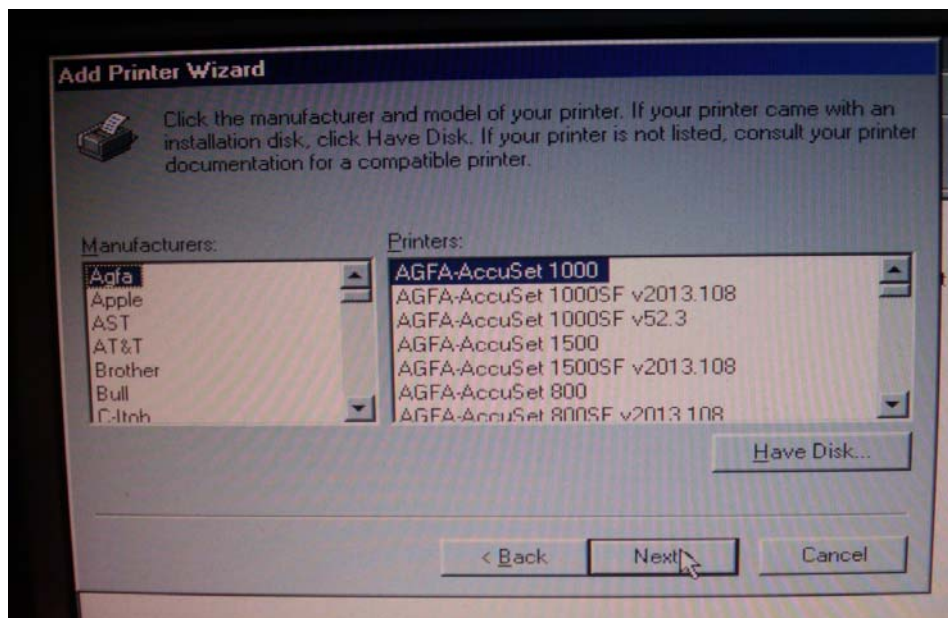


From “Printers” drop down box, select “HP LaserJet 5L PCL” and click on “Next”

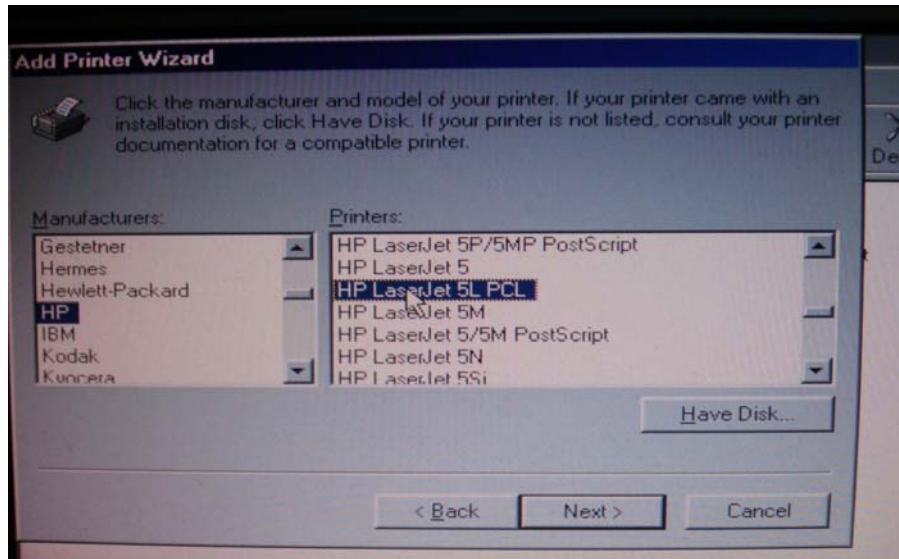


Select “LPT1: Printer Port” and click on “Next”

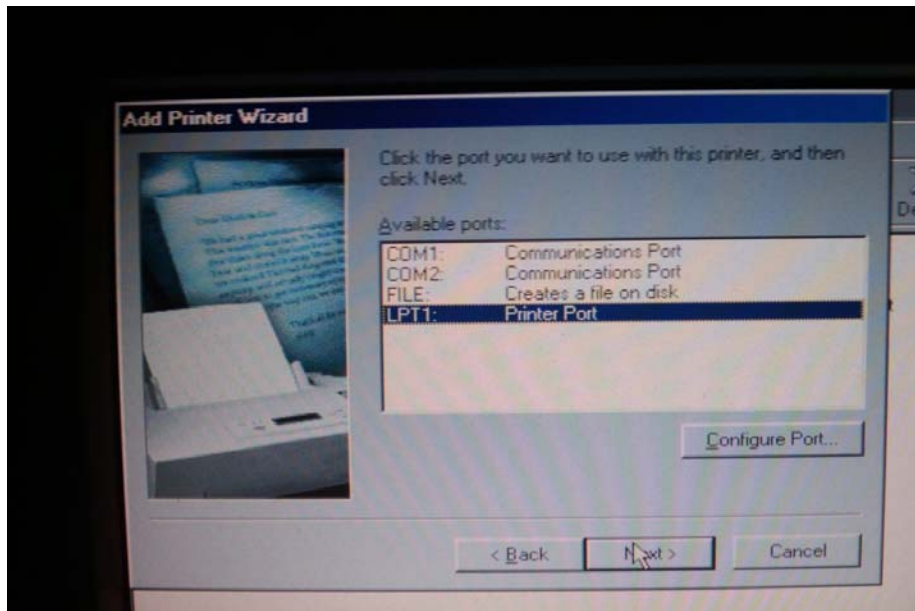


Click on the option “**Yes**” to set the new HP LaserJet as default printer and click “**Next**”

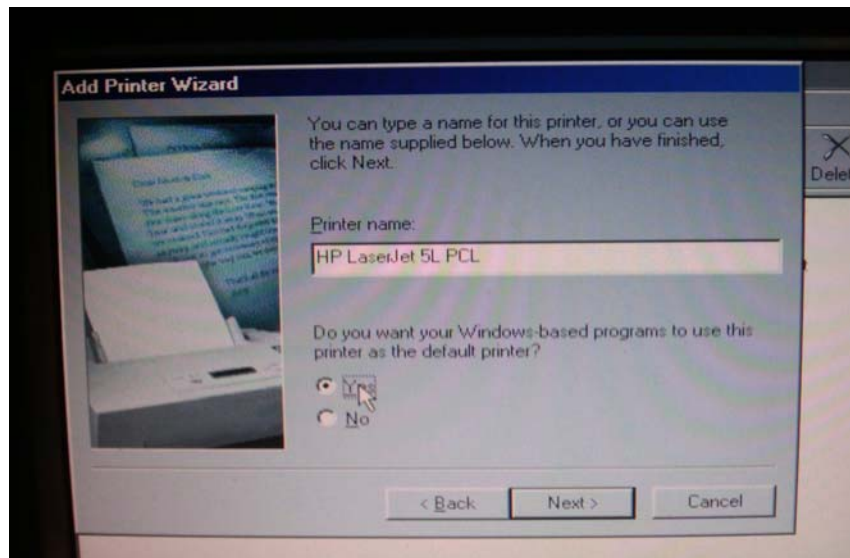


Click on the option “**No**” for disabling test page print and click on “**Finish**”

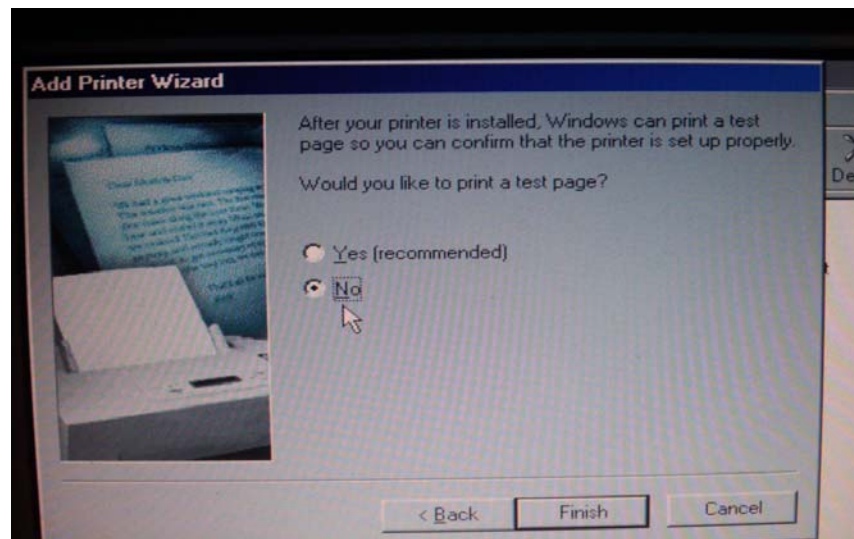


The next screen will be prompted for Windows 98 CD. Please click on **“OK”**

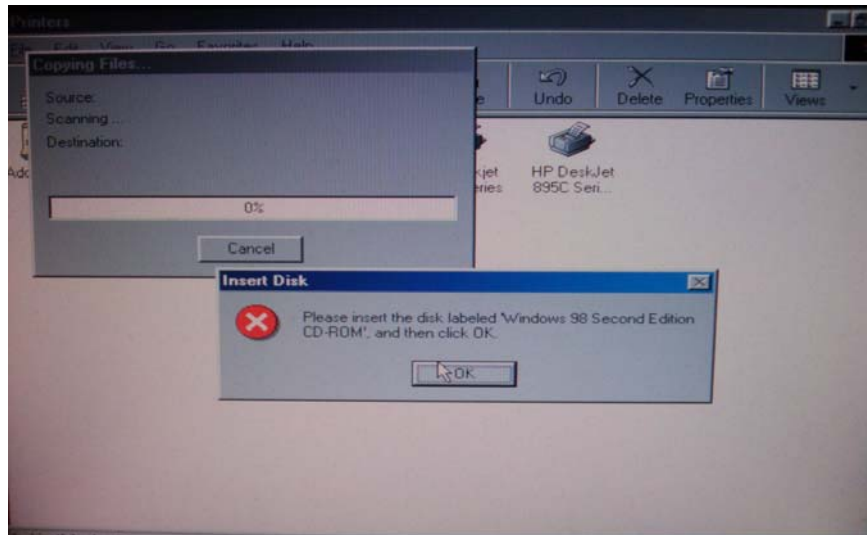


Using **“Browse”** button select files from folder **“C:\WINDOWS\options\cabs”**

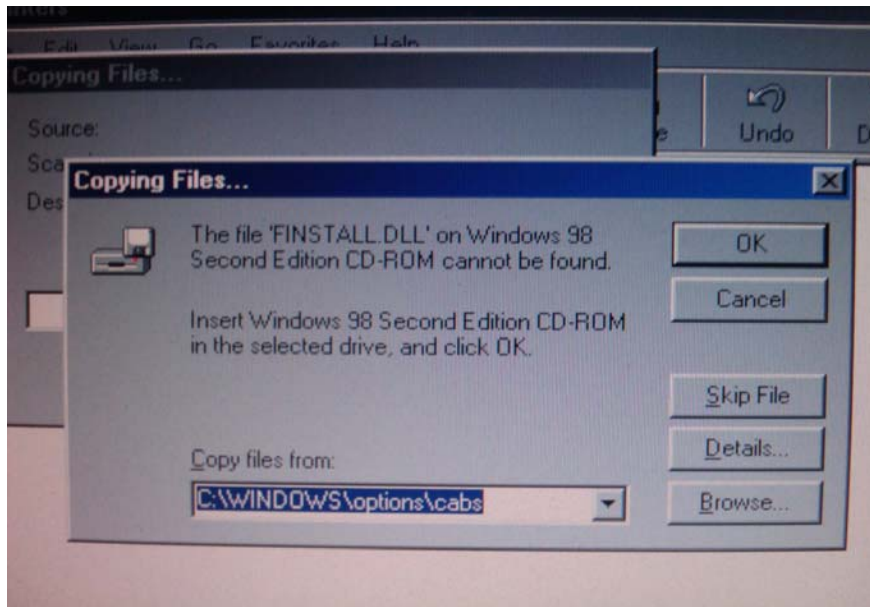

Select “**drive D**” from drives select drop down box

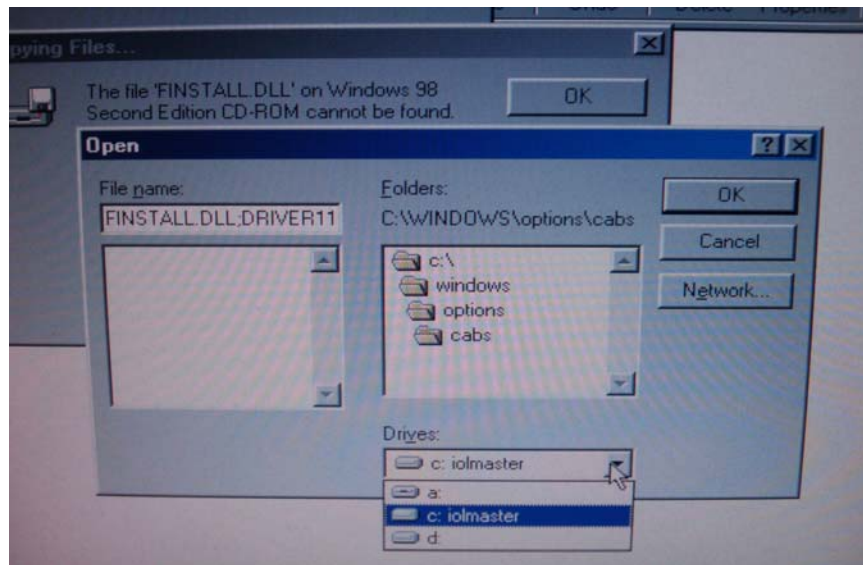


Goto folder “**SetupEng**”

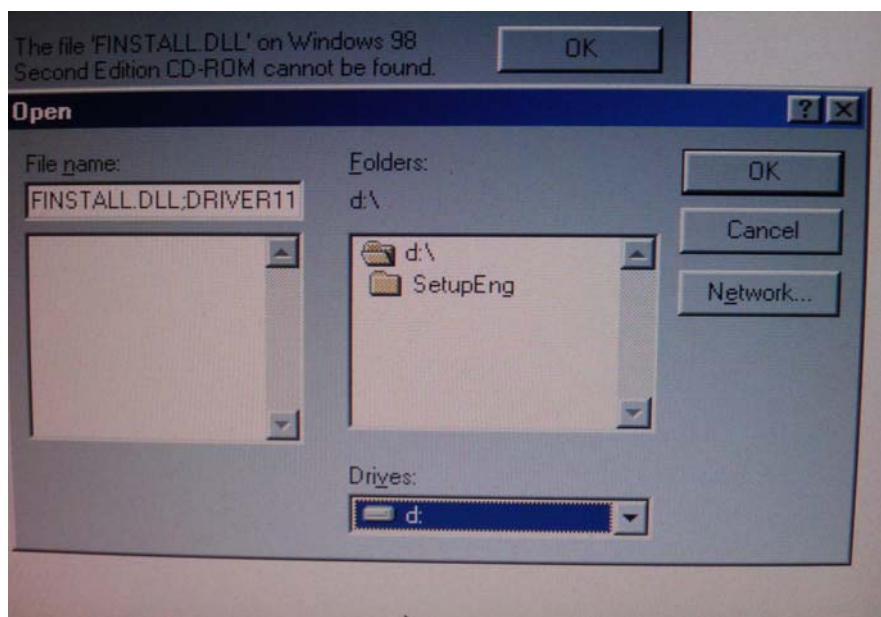





Under “SetupEng” choose “win98” folder which contains “driver11.cab” then click “OK”

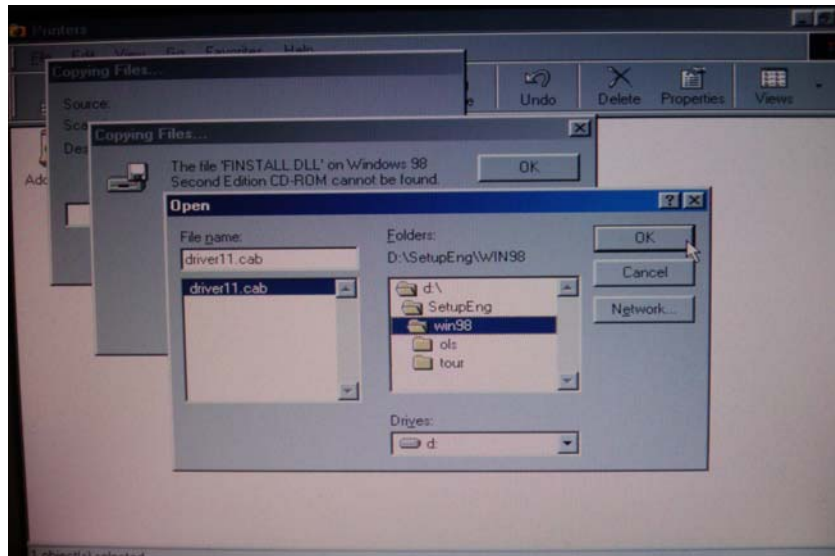


Now you can view the directory “D:\SetupEng\WIN98” as new destination and click “OK”.

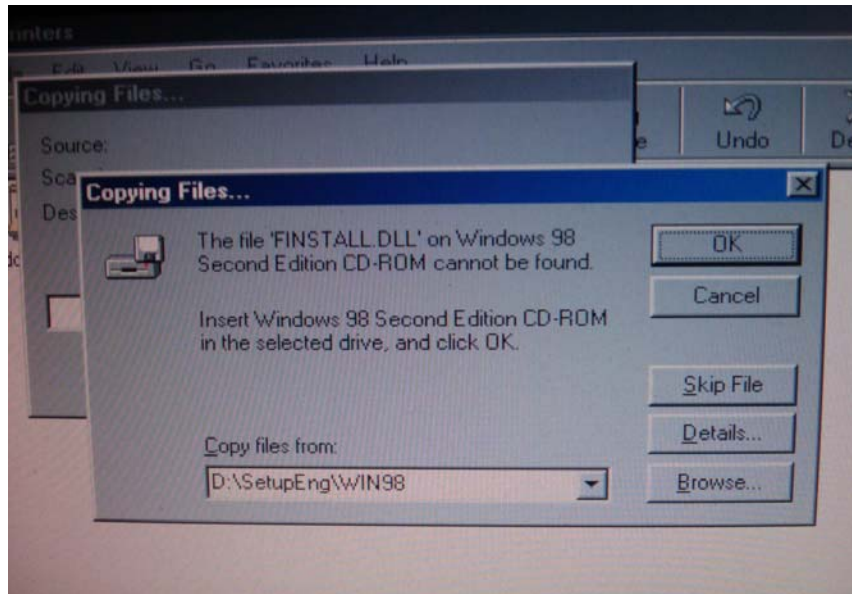


Starts copying the files from D drive to C drive.

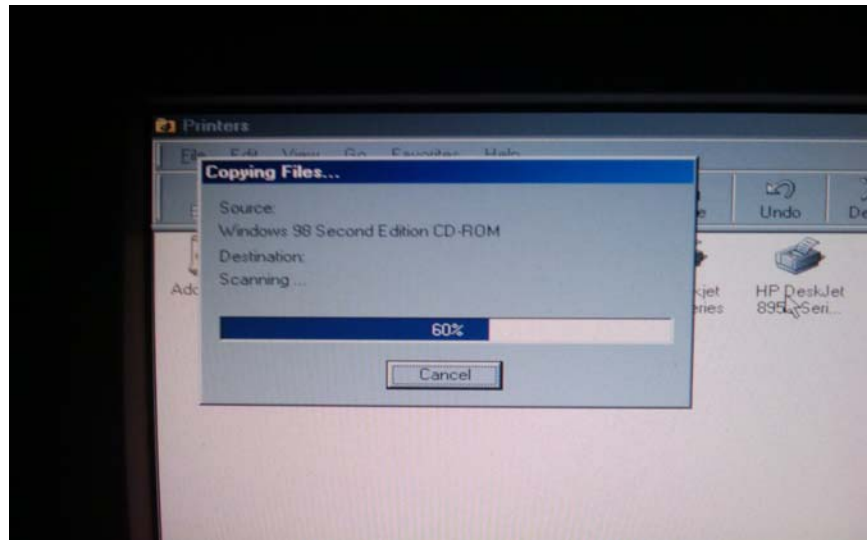

After installation of the new printer driver, you can see new printer icon as “**HP LaserJet 5L PCL**” as default printer.

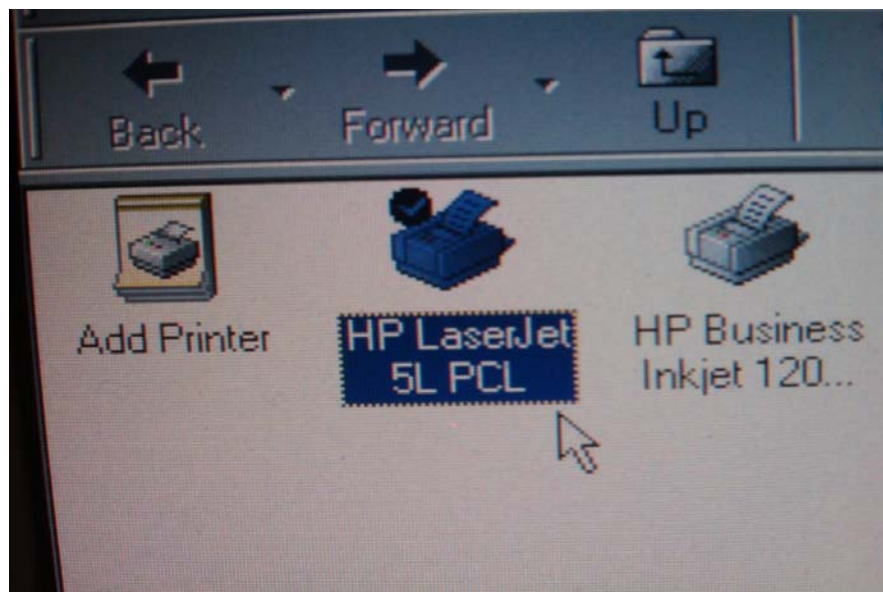

If the Print2USB box is connected, **GREEN** and **ORANGE** LED's will glow. You can choose the print option for IOL report.

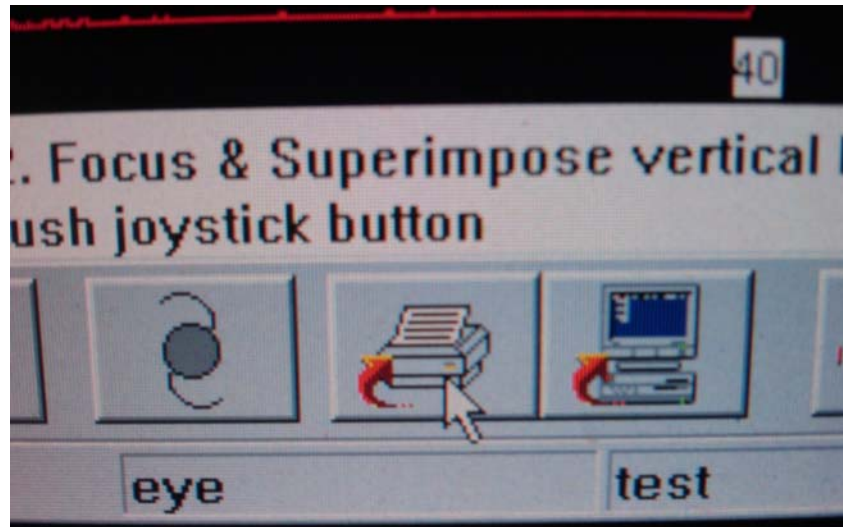

**Notes: This driver supports only monochrome print.**